

## PRESSO MASSAGE

- 🎯 Reduces water retention
- 🎯 Removes toxins of subcutaneous tissue
- 🎯 Improvement of the venous-lymphatic



## PRESSOMASS

PRESSOMASS is a presso massage device for peristaltic massage and venous-lymphatic drainage. The massage is done through a kit point with independent sectors, overlapped like a "herring bone", by means of which the pressing push gradually overlap itself without leaving interspaces during the inflating thus eliminating the danger of a back flow and improving the circulation.

This method has been designed to "drain" the limb, from the bottom to the top, slowly and continuously, respecting the venous-lymphatic system. The pressure therapy is mainly used to treat lower limbs and the areas of body that present problems of circulatory system and is recommended for all cases in which is indicated manual massage .

## A SWEET AND RELAXING TREATMENT

A cycle of sessions required 10/12 applications 2-3 times a week, with a periodic monitoring and a therapy for maintenance depending on individual needs.

The maximum pressure applied has not to exceed the minimum arterial pressure, so that the oxygenated blood can reach the tissues. We recommend not to exceed a pressure of 60 mm mercury.


## FOCUSED EFFECTS

- Improvement of the venous-lymphatic
- Reduces cellulite by increasing blood
- Removes toxins of subcutaneous tissue
- Treats the swelling of legs and ankles
- Restore a good circulation
- Reduces water retention

## THE IMPORTANCE OF "HERRINGBONE" KIT POINTS FOR LOWER LIMBS

The kit point and its structure are of fundamental importance in order to obtain optimal pneumatic compression and make sure that it cannot cause any damages. The internal structure of the kit points supplied with PRESSOMASS is called "herringbone".

For specific treatments of the venous and lymphatic system it is extremely important to use the kit points, because their pressure load overlaps smoothly without leaving interspaces where backflows can be created.

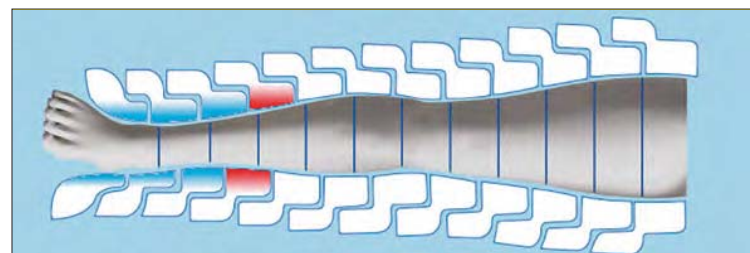

## PRESSOMASS EVO

- Trolley mounted equipment
- Color 8" touch screen display
- Membrane compressor 150 mhg
- 2 outputs
- 12 independent sectors
- Possibility to perform TOTAL BODY treatments
- Low noise membrane compresso
- Real-time reading of inflation of each sector
- Express-mode deflation: it allows a quick and silent deflation of sectors
- Warm up for all protocols: this new function introduced in software of device the treatment can be quicked, thanks to a general pre-inflation on all sectors at aprox. 30 mhg (+/-)
- Pre-set programs programs depending on the accessosies to be used

## PRESSOMASS

- Eme Estetica design
- 240x128 LCD Ghrapic display
- Trolley mounted device
- Pre-set programs
- The pressure can be setted up to 150 mm hg
- Independent sectors, each of which can be set out with differentiated work pressure
- Excellent quality of garment, with "herring bone" setting of sectors, allowing the real massage
- Smart card, where is possible to store personalized programs
- User memory
- Shelf, where is possible to store gel cosmetics
- Easy and intuitive software: possibility to load programs, possibility to personalize protocols
- Possibility to update the software



### EI10157 PRESSOMASS EVO:

Device for presso massage with 12 independent sectors, certified as CE aesthetic use

### EI101555 PRESSOMASS:

device for presso massage with 9 independent sectors, certified certified as CE aesthetic use

## PRESSOMASS 707

- Desk device
- Color 6" touch screen display
- Pre-set programs
- Independent sectors, each of which can be set out with differentiated work pressure
- Smart card, where is possible to store personalized programs
- User memory
- Easy and intuitive software: possibility to load programs, possibility to personalize protocols
- Possibility to update the software



**E110156 PRESSOMASS 707:**  
Desk device for presso massage with 7 independent sectors, certified as CE aesthetic use.

### ...SOME TIPS FOR YOU!

In order to get the best results with our equipment for press-massage PRESSOMASS, we would like to suggest its use in combination of:

COMBINATION	ULTRACAV
SEQUENCE	After the treatment with Pressomass
TARGET	To improve the lymphatic drainage of the adipocytes and of the other substances generated by their break





## Presso massage Pressomass Evo Pressomass Pressomass 707

### Features and cutting-edge technology

- Graphic display / Touch screen
  - Pre-set protocols
  - Independent sectors, each can be set with different working pressure
  - Excellent quality of garments, with "herring bone" setting of the sectors,
- allowing the real massage
- Easy and intuitive software: possibility to load programs
  - possibility to personalize protocols
  - Software Update

### Technical sheet

	PRESSOMASS EVO	PRESSOMASS	PRESSOMASS 707
Main power supply	230 Vac, 50-60 Hz + 10%, (115 Vac, 50-60 Hz+10% - upon request)	230 Vac, 50-60 Hz + 10%, (115 Vac, 50-60 Hz+10% - upon request)	230 Vac, 50-60 Hz + 10%, (115 Vac, 50-60 Hz+10% - upon request)
Display	8" touch screen color display	Grafico b/n 240x128 pixel	6" touch screen color display
Programmable treatment time	Till 99 minutes	Till 99 minutes	Till 99 minutes
Number of independent sectors	12	9	7
Maximum pressure	150 mm Hg	150 mm Hg	150 mm Hg
Stored protocols	21	21	21
Storable protocols on user's memory	200	50	100
Storable protocols on USB	4 Gb		
Storable protocols on Card		50	200
Dimensions	cm 40 x 28 x 104 H	cm 39x30x92 H	cm 39x18 x28 H
Weight	30 kg	30 kg	9 kg

### Supplied Accessories

CODE	DESCRIPTION	PRESSOMASS EVO 12 SECTORS	PRESSOMASS 9 SECTORS	PRESSOMASS 707 7 SECTORS
ACC811/E	Total Body (lower and upper limb kit point)			
ACC1317/2/E	USB key	•		
ACC810/E	12 sector lower limb kit point			
ACC536/E	9 sector lower limb kit point		•	
ACC535	5 sector upper limb kit point (gray)			
ACC535/E	5 sector upper limb kit point (bordeau)		•	
ACC534	7 sector lower limb kit point (gray)			•
ACC604	3-shelf trolley			
CONT72	Transport bag in TNT fabric			
	DVD Welcome Pack	•	•	•
	Smart Card	•	•	•